



HKN-93CS

INFANT RADIANT WARMER


Operator's manual

www.chinadavid.cn



HKN-93CS Simple Operating Instructions

1. Introduction of Controller Interface

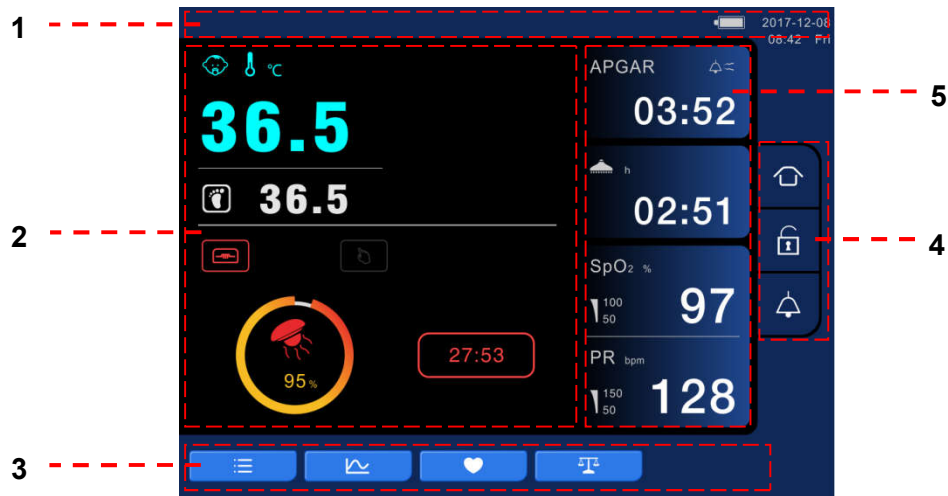

1. Status Bar

Alarm information display, Power failure alarm battery indication, RTC real time

2. Temperature control area

Skin temperature 1, Skin temperature 2, Current equipment operating mode, Heater heating power ratio

3. Menu Bar

Menu Bar, Curve function, oximetry function, Weighing function


4. Function key

Home key, Unlock/Lock key, Silence key

5. Module function

Timing function, Phototherapy function, oximetry function

2. Temperature control area

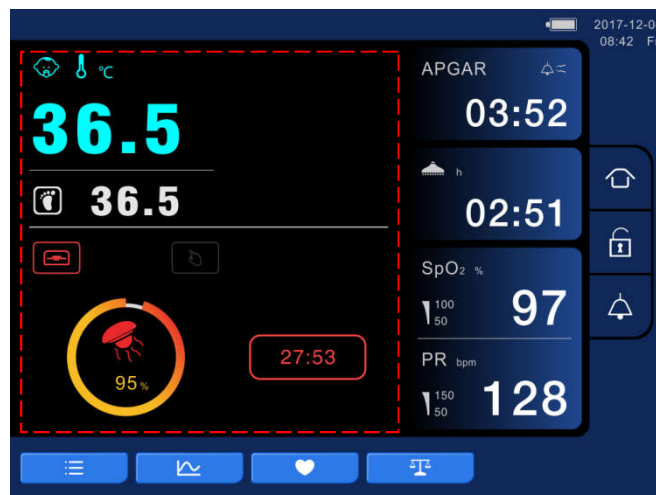

Click the red frame selection area to enter the temperature setting interface, select the operating mode of the device, and set the temperature control state.



Baby mode



Manual mode



Pre-warm mode



Increase key



Decrease key




> 37 °C setting key



OK key

3.Module function

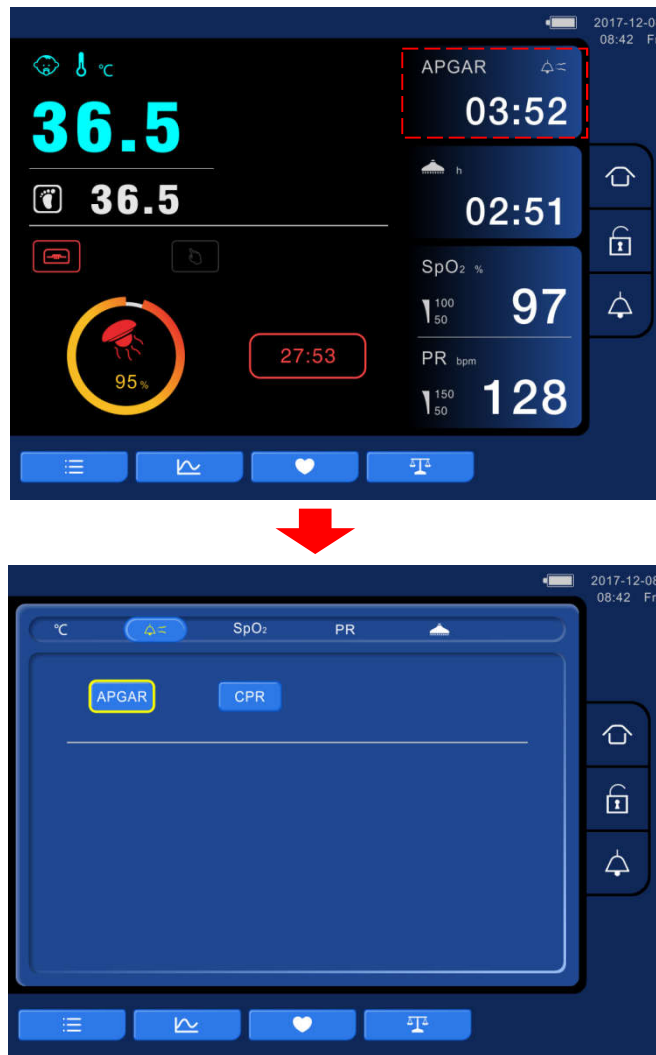


Click the red frame selection area to turn on or off the current timing mode, click temperature control area enter the timing setting interface, select the timing mode of the device.

APGAR


APGAR timing


CPR

CPR timing

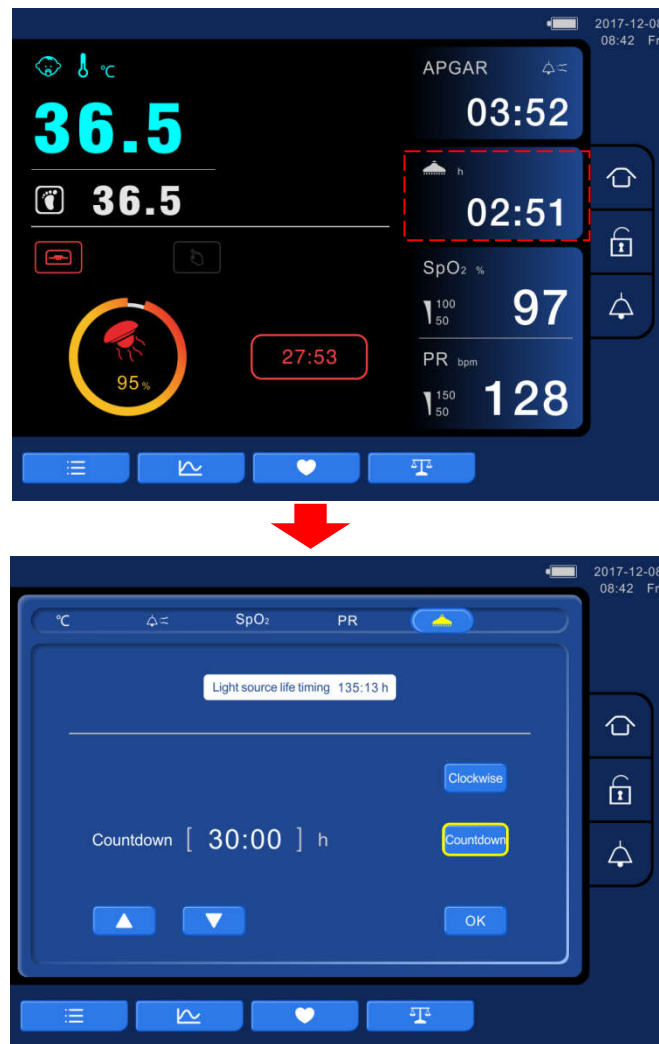



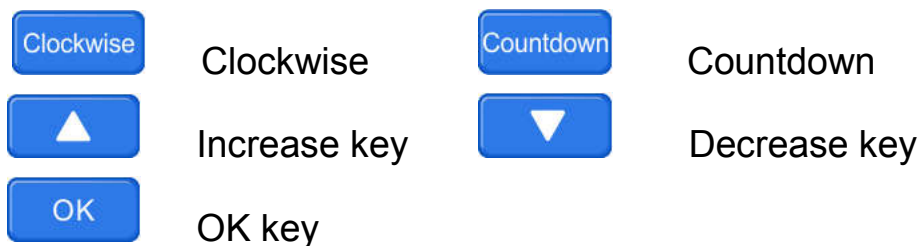

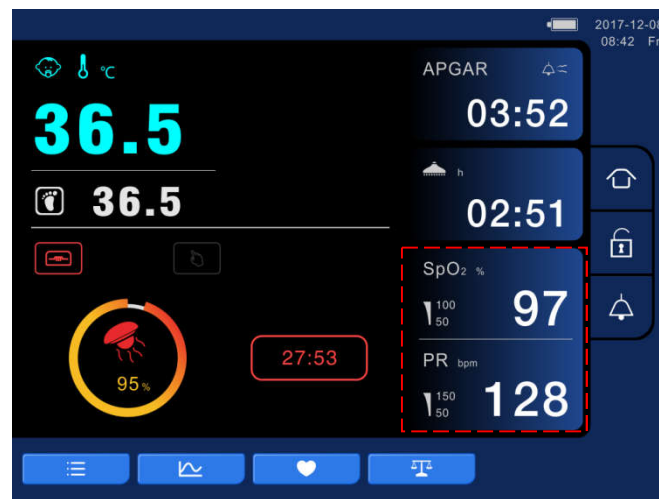


Click the red frame selection area to turn on or off the phototherap light, click temperature control area enter the phototherapy function setting interface, select the control method of the phototherapy function.





Click the red frame selection area to enter the SpO₂ setting interface, which can control the on/off of pulse oximeter monitoring module, and set the high/low alarm limit of oximetry function.


	ON		OFF
	High alarm limit key		Low alarm limit key
	Increase key		Decrease key
	OK key		

4.Menu Bar




Click the curve key to view the data tendency chart of device, including temperature data, heater heating ratio data, SpO₂ data, and pulse rate data. Click the increase/decrease key to set the data cycle of tendency chart.

	Temperature		Radiation head heating power
	Oxygen saturation		Pulse rate
	Increase key		Decrease key
	Skin temperature 1		Skin temperature 2
	Tendency chart		Tendency table

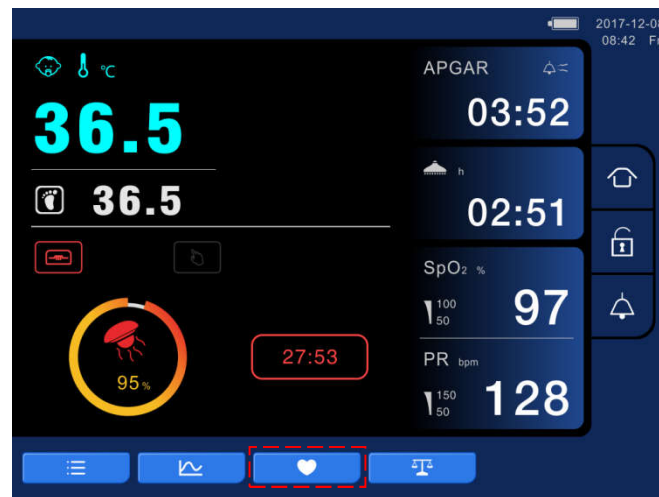

Click the oximetry key to enter the oximetry setting interface. Drag the navigation bar to set sensitivity, average time, and FastSAT .

	SENS setting		Normal mode
	APOD mode		Max mode
	AVG Time setting		FastSAT setting

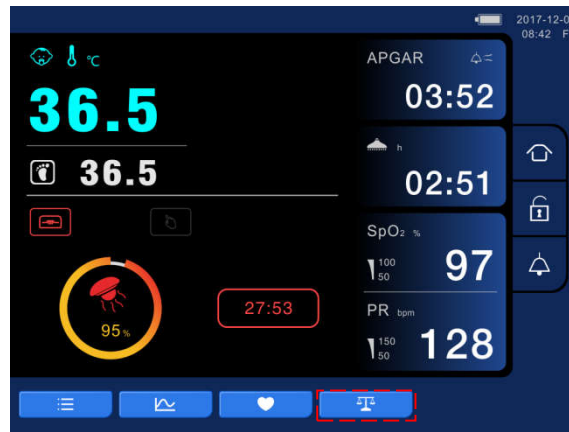

Click the baby scale key to enter the weighing interface, and drag the upper navigation bar to switch between the tendency chart and tendency table of the weight data. Click the weigh key, the baby scale function performs zero calibration, the system prompts the operator to pick up the baby 10S, then put the baby down after 10S, the weight area shows the current weight of the baby, the operator can click the record key to save the current data. If you have any doubt about the current weight, you can click the reweigh key to recalculate

the current weight; click the last weigh key to view the previous weight data, which is convenient for the operator to compare the differences between the two groups of data.



Weigh



Reweigh



Last weighing



Record



Tendency chart



Tendency table

SPECIAL STATEMENT: All of the content in the manual is checked carefully, if there is any error or content of printing misunderstanding, our company retains finally explanation of this card-usage.

NOTE: The product's appearances maybe differ from the one in this manual, but it dose not affect the capability of product. Please understand if it brings you troubles.

NINGBO DAVID MEDICAL DEVICE CO., LTD.

**ADD: NO. 2, KEYUAN ROAD, SHIPU SCIENCE AND TECHNOLOGY PARK, XIANGSHAN, 315731
NINGBO, ZHEJIANG PROVINCE, PEOPLE'S REPUBLIC OF CHINA**

FAX: +86-574-65962111

**MARKETING CENTER: BLDG NO 5. NINGBO SMART PARK, #98 CHUANGYUAN ROAD,
YINZHOU DISTRICT, NINGBO, CHINA 315048**

FAX: +86-574-87801111, 87803714

TEL: +86-574-87800008, 87800007

E-MAIL: sales@nbdavid.com

WEBSITE: <http://www.chinadavid.cn>

**EUROPEAN REPRESENTATIVE: SHANGHAI INTERNATIONAL HOLDING CORP. GMBH
(EUROPE)**

ADD: EIFFESTRASSE 80, 20537 HAMBURG, GERMANY

TEL: +49-40-2513175

FAX: +49-40-255726

EDITION/REVISION: SJKF-01

DATE OF ISSUE :2018/04